

EPSA[®] *safety first*


CYSTOMETRY CATHETER

CYSTOMETRY CATHETER

CODE: 9504 adult - 9503 pediatric

1 PVC double-lumen catheter. Filling lumen with two orifices in the distal end with blue Luer connector. Lumen for registration with one orifice at distal end with crystal colour-less Luer connector. Introduction marks.

STERILE, SINGLE PACK, SINGLE USE.

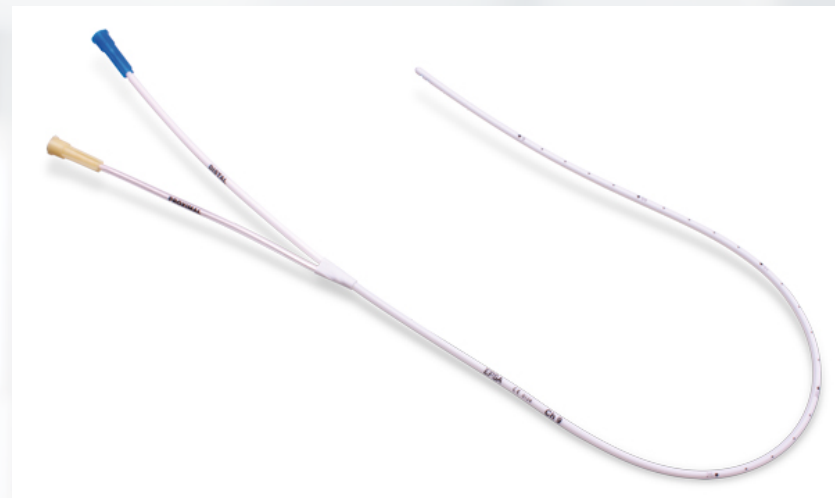

Available Sizes

Pediatric Ch 6,0 / 7,0 / Length. 32 cm.

Adult Ch 6,0 / 7,0 / 9,0 / Length. 52 cm.

Recommendation for use

Recommendation of a professional on the use of the technology.

The Cystometry Catheter is indicated for urodynamic studies of the lower urinary tract. Cystometry is a diagnostic test used to measure intravesical pressure in relation to the volume of bladder fluid.

It allows assessing the existence of voiding dysfunction in the stages of bladder filling, storage and / or emptying.

It shows the activity of the detrusor muscle and the sensitivity and capacity for contraction and relaxation of the bladder wall at rest and in relation to some stimuli.

