



EPSA[®] *safety first*


UMBILICAL CATHETER

UMBILICAL CATHETER

CODE: 0631

1 PUR catheter, graduated in cm., continuous longitudinal radiopaque line. Polished tip, atraumatic. Luer connector at the proximal end.
1 polyethylene protective packaging.

STERILE, SINGLE PACK, SINGLE USE.

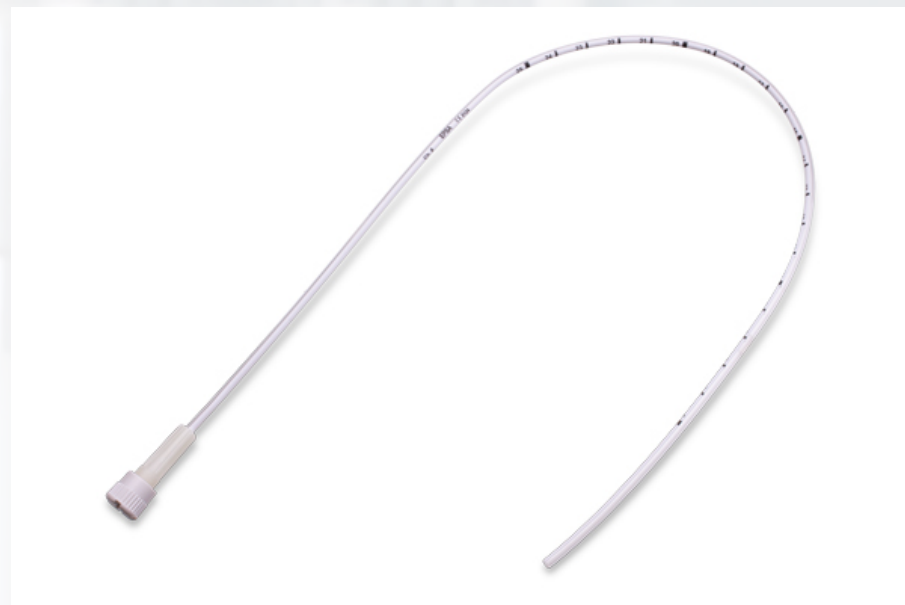

Available Sizes

Ch 2,5 / 3,0 / 3,5 / 4,0 / 5,0 / 6,0 / 7,0 / 8,0.

Long. 40 cm

Recommendation for use

Recommendation of a professional on the use of the technology.

The Umbilical Catheter is indicated for catheterization of the vessels (arteries or veins) of the newborn's umbilical cord. It is used in numerous pathologies of the neonate.

– Umbilical arterial catheterization (UAC): the purpose of this catheter is to continuously monitor blood pressure and allow frequent blood sampling for analysis without having to perform repetitive punctures. It is used most often if the newborn has severe lung disease that requires mechanical ventilation or blood pressure disturbances that require the use of special medication.

– Umbilical venous catheterization (UVC): the purpose of this catheter is to allow the administration of fluids and medications without having to replace the IV repeatedly. It is used most often in immature newborns, with intestinal problems that prevent adequate nutrition and those who require the provision of special medication.

